


Mectalif ANTERIOR

ANTERIOR LUMBAR INTERBODY FUSION DEVICE

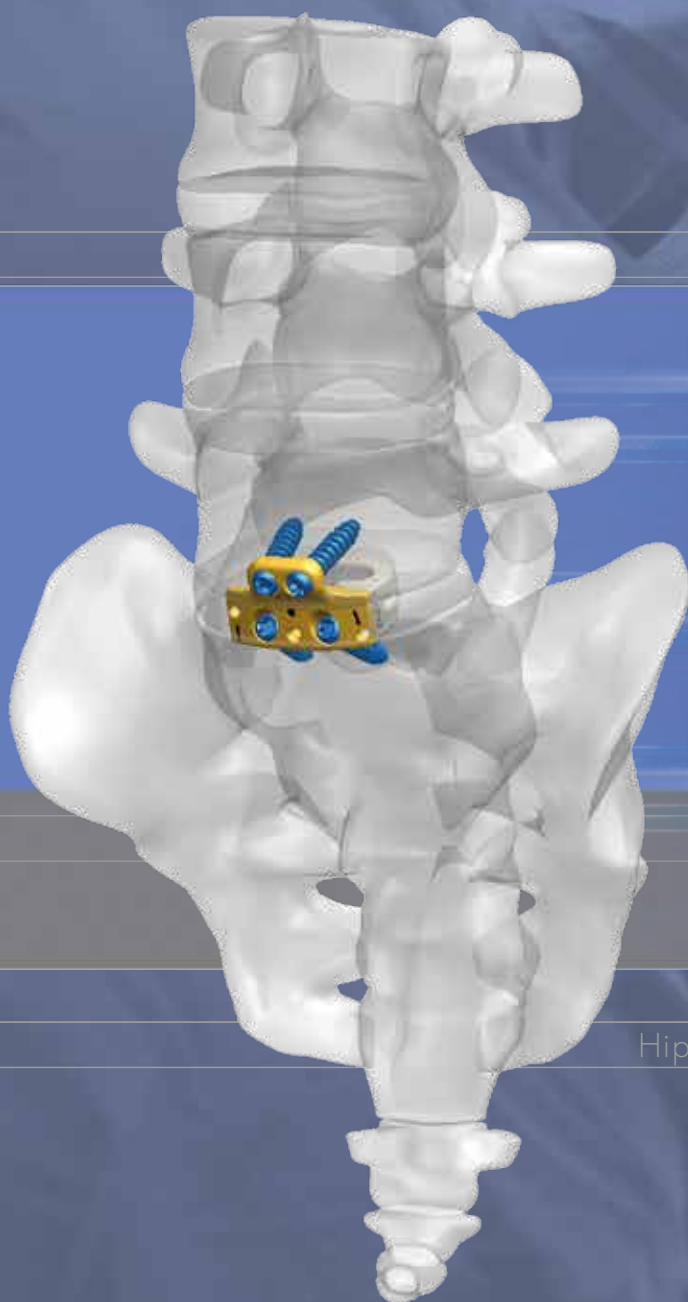

Brochure

Hip

Knee

Spine

Navigation

MECTALIF ANTERIOR HYBRID

The MectaLIF **Hybrid** Interbody Fusion Device merges the mechanical features of the **Flush** zero anterior profile plate with the **Long** enhanced stability construct, thus providing added value in Anterior lumbar fixation.



MECTALIF ANTERIOR HYBRID RATIONALE


The MectaLIF **Hybrid** Interbody Fusion Device provides an aided solution for L5-S1 implantation where additional stability is required and low profile construct is necessary to deal with difficult tilted anatomies.



The **Flush** caudal profile avoids impingement and irritation of the soft tissues.



The **Long** cranial profile provides high level stability in extension and torsion.

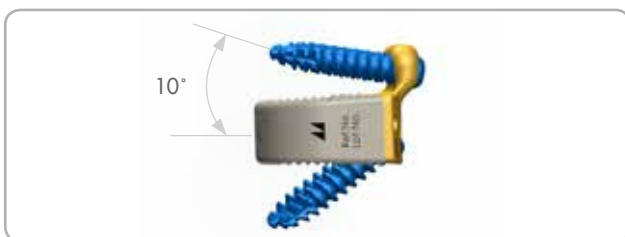

The low-angle 10° screws offer easier insertion in challenging anatomies avoiding interference with the pubic symphysis. Low angle trajectory allows for screw insertion with straight instruments and a simple surgical technique.



Diverging screws provide enhanced stability and increased pull-out strength.

REFERENCES	
Hybrid Plate Height	Ref. code
H.12 mm	03.30.232
H.14 mm	03.30.233
H.16 mm	03.30.234
H.18 mm	03.30.235