


# **Mectalif**<sup>®</sup> ANTERIOR

ANTERIOR LUMBAR INTERBODY FUSION DEVICE

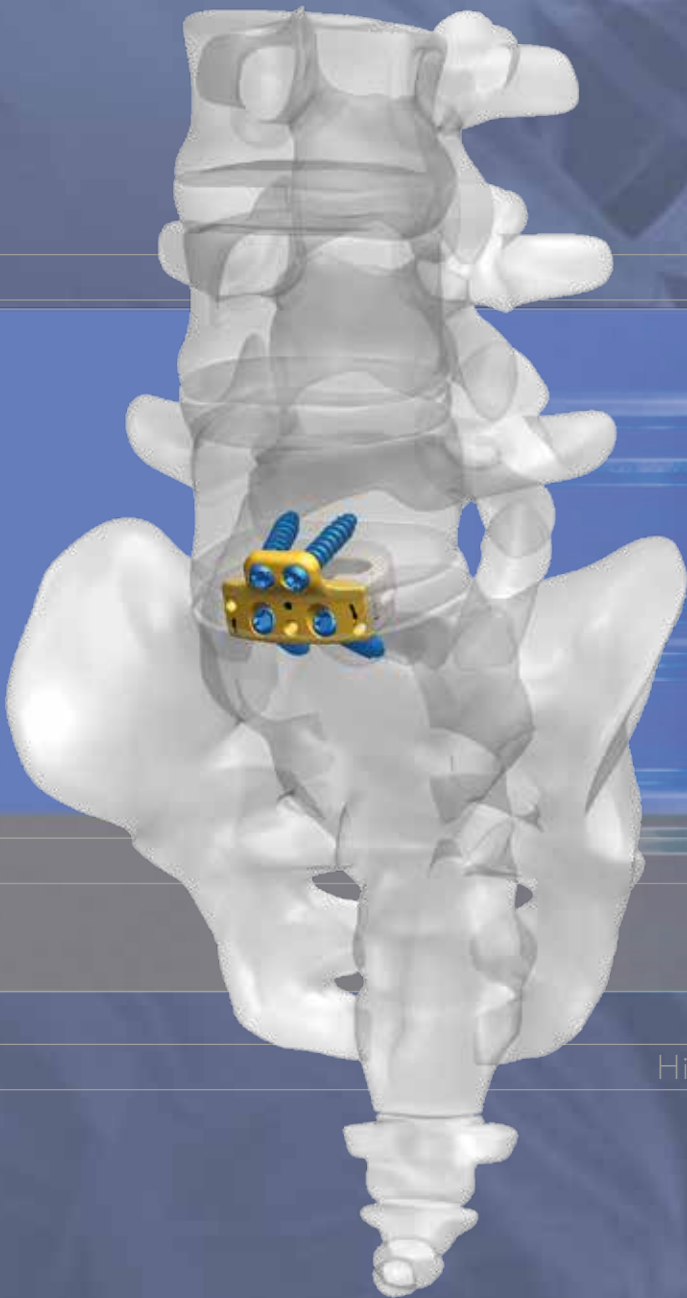

Brochure

Hip

Knee

Spine

Navigation

## MECTALIF ANTERIOR HYBRID

The MectalIF **Hybrid** Interbody Fusion Device merges the mechanical features of the **Flush** zero anterior profile plate with the **Long** enhanced stability construct, thus providing added value in Anterior lumbar fixation.



## MECTALIF ANTERIOR HYBRID RATIONALE


The MectalIF **Hybrid** Interbody Fusion Device provides an aided solution for L5-S1 implantation where additional stability is required and low profile construct is necessary to deal with difficult tilted anatomies.



The **Flush** caudal profile avoids impingement and irritation of the soft tissues.



The **Long** cranial profile provides high level stability in extension and torsion.

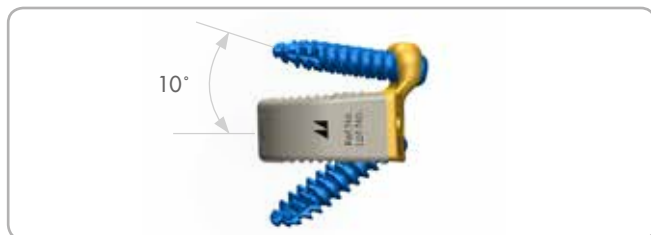

The low-angle 10° screws offer easier insertion in challenging anatomies avoiding interference with the pubic symphysis. Low angle trajectory allows for screw insertion with straight instruments and a simple surgical technique.



Diverging screws provide enhanced stability and increased pull-out strength.

REFERENCES	
Hybrid Plate Height	Ref. code
H.12 mm	03.30.232
H.14 mm	03.30.233
H.16 mm	03.30.234
H.18 mm	03.30.235

All trademarks and registered trademarks are the property of their respective owners.

