


MectaLock PEEK RC
J-LOCK PEEK SUTURE ANCHOR

MectaLock C RC
J-LOCK COMPOSITE SUTURE ANCHOR

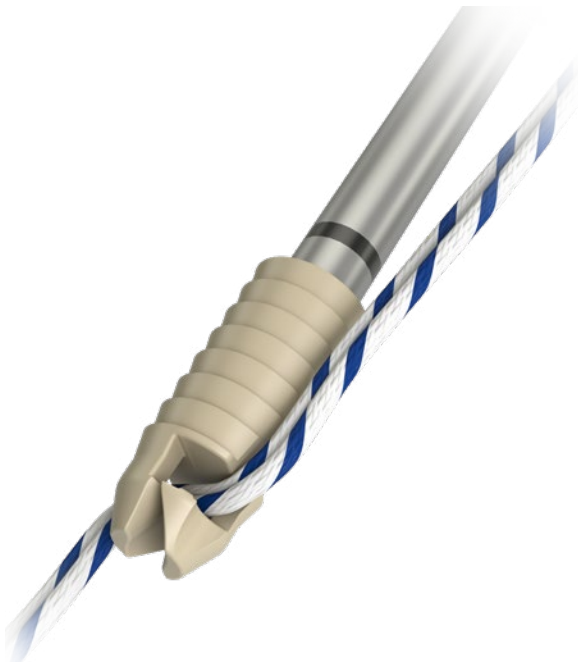

KNOTLESS ANCHOR FOR ROTATOR CUFF REPAIR

NO DRILLING REQUIRED

Pilot hole is created with a punch
Anchor is malleted in afterwards

TECHNICAL DATA

IMPLANT LENGTH: 15mm
Available in 3 sizes:
4.5mm/ 5.5mm/ 6.5mm

SUITABLE FOR LATERAL ROW ROTATOR CUFF REPAIR

SUTURE EYELET INTEGRATED IN ANCHOR

MATERIAL

2 different materials available:

- MectaLock PEEK RC


- MectaLock C RC
Composite, Osteoconductive

70%	30%
PLDLLA	+ β -TCP



Reusable punch to create pilot hole

ANCHOR BODY
The interference between the surrounding bone and the ribbed anchor body increases pull-out strength

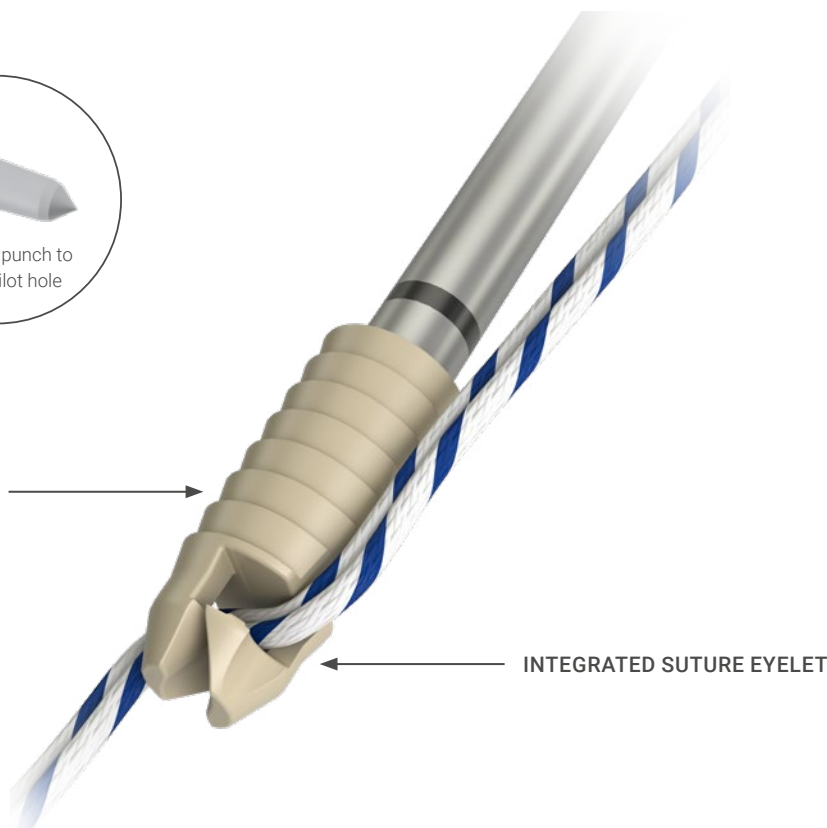

INTEGRATED SUTURE EYELET

ORDERING INFORMATION

MECTALOCK PEEK RC			
Order Number	Description	Order Number	Description
05.10.004	MectaLock PEEK - Ø4.5 - Short driver	05.10.006	MectaLock PEEK - Ø6.5 - Short driver
05.10.005	MectaLock PEEK - Ø5.5 - Short driver		

MECTALOCK C RC			
Order Number	Description	Order Number	Description
05.10.104	MectaLock C - Ø4.5 - Short driver	05.10.106	MectaLock C - Ø6.5 - Short driver
05.10.105	MectaLock C - Ø5.5 - Short driver		

DEDICATED INSTRUMENTS			
Order Number	Description	Order Number	Description
05.10.10.0021	Punch - Ø4.5	05.10.10.0023	Punch - Ø6.5
05.10.10.0022	Punch - Ø5.5		

RELATED PRODUCTS

Order Number	Description	Order Number	Description
05.15.10.0015	FastShuttle - Suture Passer & Retriever	05.14.10.0001	Shoulder Cannula System Ø8.5x70mm
05.15.10.0016	FastShuttle - Suture Passer	05.14.10.0002	Shoulder Cannula System Ø8.5x90mm
05.15.10.0017	FastShuttle - Disposable Needle - STERILE (5ea)	05.14.10.0035	Shoulder Cannula System Ø8.5x70mm for Ø4mm Switching Stick
05.10S.001	Shoulder General Set	05.14.10.0036	Shoulder Cannula System Ø8.5x90mm for Ø4mm Switching Stick



[SPORTSMED.MEDACTA.COM](https://www.sportsmed.medacta.com)

Distributed by **Medacta International SA** | Strada Regina - 6874 Castel San Pietro - Switzerland | Phone +41 91 696 60 60 | Fax +41 91 696 60 66 | info@medacta.ch

Find your local dealer at: medacta.com/locations. All trademarks and registered trademarks are the property of their respective owners. This document is not intended for the US market. The devices described in this document are not available in all the Countries in which Medacta operates. Please verify approval of the devices described in this document with your local Medacta representative.

ref: 99.125SMS.1SALES-S rev. 00 Last update: April 2022